



**Battery Power  
Solutions**

Batteries to Power Your World!

**Freedom<sup>TM</sup>V2**  
CPAP BATTERY

# Auto Charge DC Cable



ITEM # DC.CABLE-V2



# Auto Charge DC Cable Instructions

Charge the Freedom V<sup>2</sup> CPAP Battery via the power outlet of a vehicle. The vehicle must be running during charging to ensure the battery receives the proper voltage to successfully charge to 100%.

## Specifications


Input/Output: 12V DC | Charge Time: 3-4 hrs | Length: 53" | Weight: 3 oz

## Single Battery Charging

- Plug the Auto Charge DC Cable into the DC IN port of the battery.
- Turn the vehicle on and plug the Auto Charge DC Cable into the power outlet of the vehicle to begin the charging process.

## Dual Battery Charging

- Plug the Bridge Cable into the DC out port of battery A and then into the DC IN port of battery B. **Battery A must be charged to at least 15% but battery B can be at any charge level including 0%.**
- Plug the Auto Charge DC Cable into the DC IN port of battery A.
- Turn the vehicle on and plug the Auto Charge DC Cable into the power outlet of the vehicle to begin the charging process.

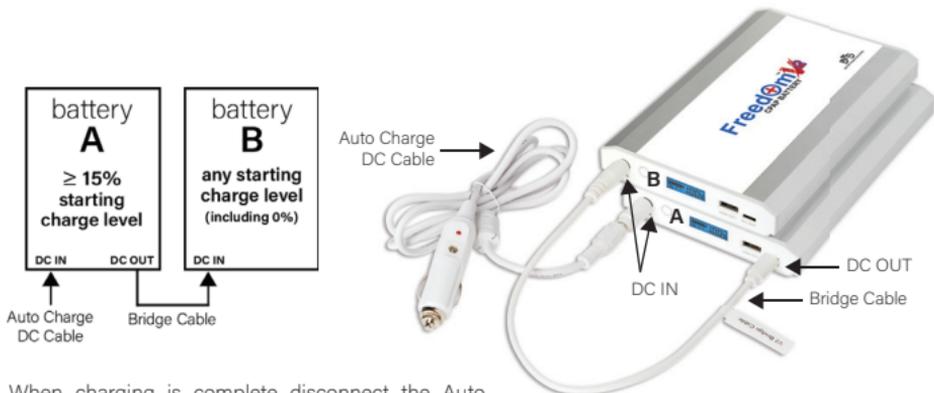

When charging is complete disconnect the Auto Charge DC Cable from the vehicle and the battery.